ENGLISH4

ENGLISH

FRANÇAIS 10

FRANÇAIS

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

Foot SPA Massager FBM-311

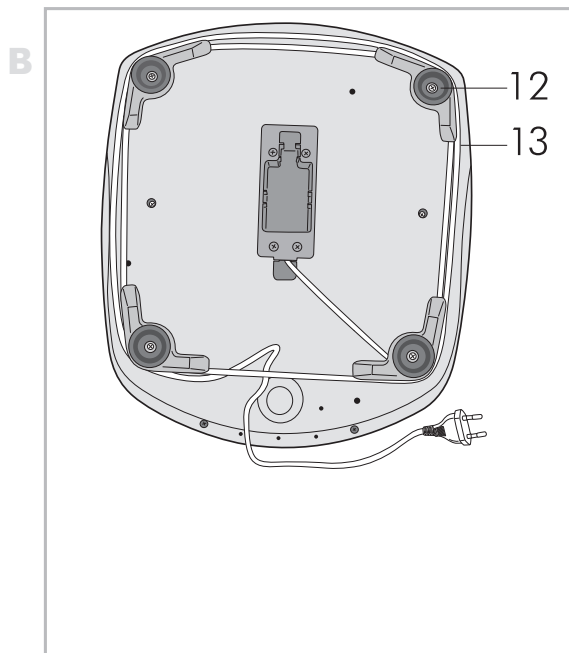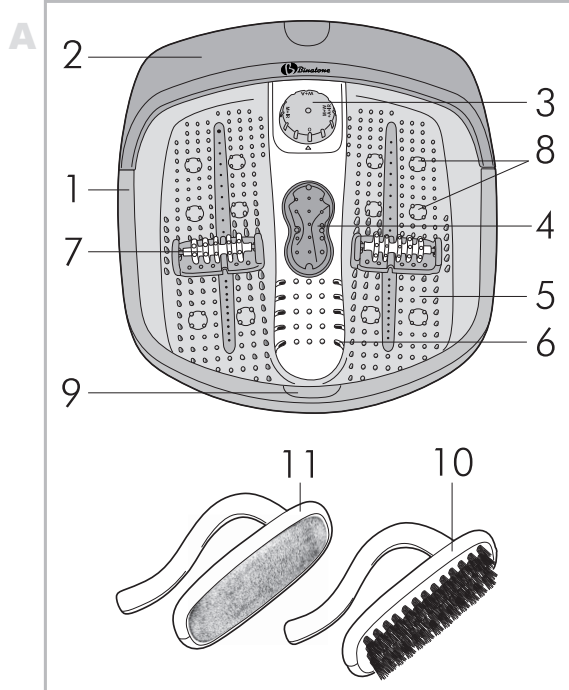Foot SPA Massager FBM-311

Please read the safety instructions carefully before using this appliance. **Please retain the instructions for future reference.**

SAFETY MEASURES

When using your appliance, the following basic safety precautions should always be followed:

- ▶ Before using for the first time, check that the appliance voltage is the same as your local supply.
- ▶ Only use the appliance for the purposes stipulated by the instructions.
- ▶ Only use attachments supplied with the equipment.
- ▶ To avoid electric shocks, **do not immerse the appliance, plug or power cord in water or any other liquids.** Do not place the appliance where it could come into contact with water or other liquids. If the appliance comes in contact with water, immediately disconnect it from the mains. Do not use the machine again until it has been serviced by **an authorized service center.**
- ▶ Do not let children use the appliance without supervision.
- ▶ **If the power cord plug does not fit to your wall socket, consult a qualified electrician.**
- ▶ **Use of various adapters may result in damage to the appliance and termination of warranty.**
- ▶ Always supervise the equipment while in use.
- ▶ To avoid a circuit overload, do not operate another high voltage equipment on the same circuit.
- ▶ **Always unplug before moving the appliance from one place to another, cleaning and storing.**
- ▶ Never pull the cord to disconnect the appliance, instead grasp the plug and pull it to disconnect.
- ▶ Do not allow the cord or the machine itself to come in contact with hot surfaces.
- ▶ Position the equipment and the cord so that nobody can accidentally overturn the equipment.
- ▶ Do not use the appliance if it has been dropped, or it has a damaged power cord or plug. Do not modify or repair the machine yourself.
- ▶ Do not dismantle the appliance as improper assembly might result in an electric shock during the subsequent use of the equipment. **Repairs should only be done by authorized service centers.**
- ▶ Do not use the appliance where the air may contain volatile flammable substances.
- ▶ This equipment is intended for household use only, not for commercial or industrial use.
- ▶ Please follow all the requirements of this guide.
- ▶ Do not use the appliance if it has mechanical defects (dents, cracks etc.), get it checked by the nearest authorized service center.

SPECIAL PRECAUTIONS FOR USE

ENGLISH

- ▶ Use pure water and bath salts (oils) that do not produce foam.
- ▶ Dissolve salts and oils separately in order not to cork up air feeding holes.
- ▶ Never carry water in the apparatus above the "MAX" mark.
- ▶ Do not move the apparatus or unplug from the socket during operation of the apparatus.
- ▶ Never get into the apparatus. Use it only in sitting position.
- ▶ Do not use hydromassage function without water.
- ▶ Dispose of the water through the special lip.
- ▶ To avoid fire or electric shock, do not remove leads designed for maintenance of the electrical appliance. The electrical appliance has no elements that require user service.
- ▶ Avoid ingress of water on the regulator of operation mode.
- ▶ Do not use for more than 20 minutes. After 20 minutes of continuous work, the apparatus should be switched off and allowed to cool to room temperature.
- ▶ It is not advisable to use massage apparatus if you have phlebitis, varicose veins, severe leg edema, burns or open wounds, diseases of the nervous system (for example sciatica) or during pregnancy. In the above cases, we recommend you to consult your doctor before use of the massage apparatus.
- ▶ When in use the appliance has a hot surface. Persons insensitive to heat must be careful when using this appliance.

DESCRIPTION OF APPARATUS COMPONENTS

Figure A:

1. Case
2. Cap
3. Switch for operation mode:
 - 0 – off
 - 1 – vibration, massage with infrared radiation
 - 2 – hydro-massage with heating
 - 3 – hydro-massage with vibration, heating and infrared radiation
4. Infrared lamp
5. Surface for acupunctural massage
6. Holes for air bubbles
7. Removable massage rollers for intense foot massage
8. Built-in magnets
9. Lip for water disposal
10. Brush for pedicure
11. Pumice stone

- 12. Rubber feet (**Figure B**)
- 13. Cord storage (**Figure B**)

BEFORE USE

Remove all package materials from the apparatus and clean it as indicated in this instruction.

USER INSTRUCTION

Use of the apparatus

1. Install the apparatus on a flat horizontal surface in front of a comfortable seat (chair).
2. Fill the massage apparatus with water of required temperature to the maximum level. To enhance the effect of the massage, you may add salt or bath oil to water.
3. Connect the apparatus to electric mains.
4. Sit down and put down feet in bath. **Always put in and take out feet from the bath in whilst in sitting position.**
5. Set the switch of operation mode (**3**) in required position.
6. For more intense foot massage, install massage rollers (**7**) on the surface for acupunctural massage (**5**), as shown in **figure A**.
7. After massage procedures, set the switch (**3**) in position «0» and follow instructions in the section «Maintenance and cleaning».

Description of operation modes of the apparatus

Heating

Massage apparatus heats water to 40°C, irrespective of its starting temperature. Hence, you can fill with cold water. Remember, that heating the water will take time.

Infrared radiation

Infrared radiation improves circulation in legs and relieves muscular stress. Infrared radiation leads to fast cell renewal. It has an anti-inflammatory and analgesic effect.

Magnetic field

Magnetic field provides maximum relax and stimulates circulation in legs.

Vibration massage with infrared radiation

Vibration massage relieves muscular stress and improves circulation. Vibration massage can be used both with and without filling the apparatus with water. Combination of infrared radiation and massage gives a very good therapeutic effect.

Note: time of vibration massage without water should not exceed 10 minutes, with water — 20 minutes.

Hydromassage with heating

Warm water in combination with bubble effect produces a sensation of light touch which gives a soothing effect on the nervous system as a whole. Air, fed through small holes, produces bubbles, which gently massage feet and skin, providing pleasure and relaxation.

Hydromassage with vibration, heating, and infrared radiation

This mode provides intense massage. Relieves leg fatigue and improves health state.

MAINTENANCE AND CLEANING

1. Firstly disconnect the apparatus and pull out plug from socket.
2. Then pour away the used water from massage apparatus through the lip (9).
3. Carefully wash the internal surface of the apparatus with water, using a mild cleanser and dry it, particularly after using oils.
4. Wash detachable nozzles (7) and accessories (10 and 11) with soapy warm water.
5. After use wrap cord around apparatus legs, as shown in figure B.

USEFUL ADVICE

Adding aromatherapy treatment oils enhances the hydromassage efficiency.

Problem	Oil dosage
Foot fatigue	<ul style="list-style-type: none"> • 5 drops of rosemary oil, 2 drops of lemon oil • 5 drops of Japanese mint oil • 3 drops of geranium oil, 3 drops of lemongrass oil
Insufficient blood supply	<ul style="list-style-type: none"> • 2 drops of geranium oil, 2 drops of cypress oil, 5 drops of rosemary oil
Foot sweating	<ul style="list-style-type: none"> • 8 drops of sage oil • 6 drops of lavender oil

Please see the chart below for the recommendations on how to use essential oils for alleviating foot problems.

DISPOSAL

This appliance should be disposed of in accordance with waste disposal regulations in your country and with the least damage to the environment.

Only ecologically safe materials are used in the packaging of our products and the cardboard and paper can be disposed of with normal paper waste.

AFTER-SALES SERVICE

If your appliance does not operate at all or operates poorly, call an authorized service center for help or repairs.

The warranty does not cover the following:

- ▶ Non-observance of the above safety measures and operation instructions.
- ▶ Commercial use (use outside the scope of personal household needs).
- ▶ Plugging into a mains socket with a voltage different from the voltage on the appliance rating plate.
- ▶ Improper care (ingress of liquid, dust, insects, and other foreign items and substances inside the appliance).
- ▶ Use of excessive force during appliance operation or other damages (mechanical damages resulting from incorrect use, shocks or dropping of the appliance).
- ▶ Unauthorized repairs or replacement of the appliance parts performed by an unauthorized person or service center.
- ▶ Accidents, such as lightning strike, flood, fire and other acts of God.

ENGLISH

Features may be changed by Binatone without any notice.
Warranty details are set out in the warranty card provided with the product.

Exploitation period: 3 years after the day of purchase.

Manufacture: Binatone Industries Ltd, Great Britain, Vitabiotics House, 1 Apsley Way, Staples Corner, London, NW2 7HF, UK.

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