SCARLETT™ *Vita*SPA

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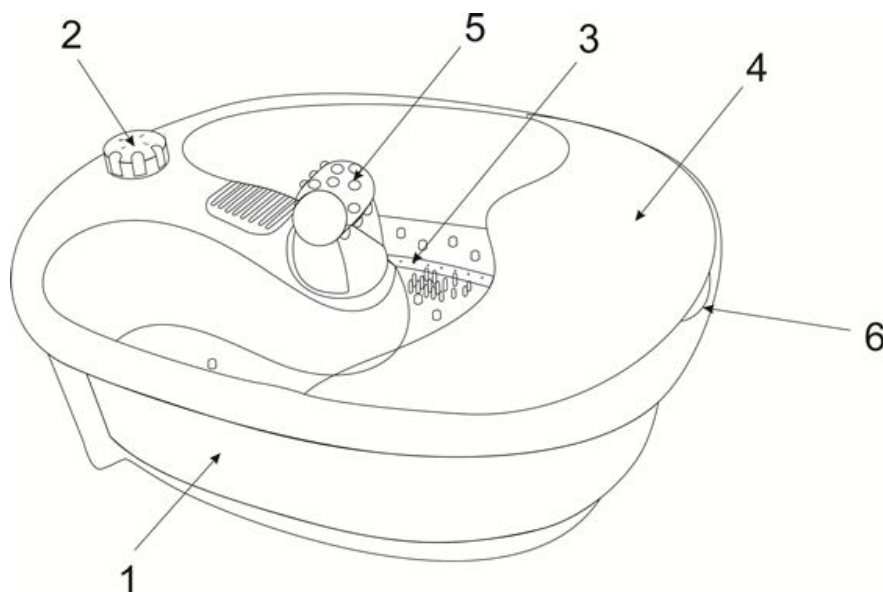
GB INSTRUCTION MANUAL

IMPORTANT SAFEGUARDS

- Please read this instruction manual carefully before use and keep in a safe place for future reference.
- Before the first connecting of the appliance check that voltage indicated on the rating label corresponds to the mains voltage in your home.
- For home use only. Do not use for industrial purposes. Do not use the appliance for any other purposes than described in this instruction manual.
- Do not use outdoors.
- Always unplug the appliance from the power supply before cleaning and when not in use.
- To prevent risk of electric shock and fire, do not immerse the appliance in water or any other liquids. If it has happened DO NOT TOUCH the appliance, unplug it immediately and check in a service center.
- Do not take this product to a bathroom or use near water.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliances by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Do not leave the appliance switched on when not in use.
- Do not use other attachments than those supplied.
- Do not operate after malfunction or cord damage.
- Do not attempt to repair, adjust or replace parts in the appliance. Check and repair the malfunctioning appliance in the nearest service center only.
- Keep the cord away from sharp edges and hot surfaces.
- Do not place or store the appliance where it can fall or be pulled into a tub or sink.
- Do not pull, twist, or wrap the power cord around the appliance.
- The massager is designed for intermittent use: it can be applied without interruption not longer than **10** minutes, then the motor must be allowed to cool.
- Do not use the massager if you are pregnant, or suffer from thrombosis, open wounds, broken skin and varicose veins. Consult your doctor first.
- Do not attempt to plug in or unplug the unit while feet are in water.
- Do not stand up or place any objects in the unit, use it only when you are sitting.

- Immediately interrupt procedure if you feel discomfort, or if any pain or irritation results; never use for swollen or inflamed feet, or if feet have any skin eruptions. First, obligatory consult your physician.
- If a unit is leaking, the device can not be used.
- **DO NOT PLUG THE DEVICE IN THE SOCKET IF YOUR FEET ARE INSIDE THE MASSAGE TUB.**
- Do not attempt to repair, adjust or replace parts in the appliance. Repair the malfunctioning appliance in the nearest service center.
- If the product has been exposed to temperatures below 0°C for some time it should be kept at room temperature for at least 2 hours before turning it on.
- The manufacturer reserves the right to introduce minor changes into the product design without prior notice, unless such changes influence significantly the product safety, performance, and functions.

GENERAL INFORMATION

- Regular foot massage dramatically improves blood circulation, relieves muscle tension, increases vitality, allows making the skin of your feet softer.
- Built-in heating element maintains the temperature throughout the whole massage session.
- 36 holes for air in the massage surface provide a relaxing bubble bath.
- The bottom of the bath is ergonomically designed and covered with massage protuberances that are needed for acupressure. The massage strength can be easily changed by pressing the feet against the massage surface or by easing the pressure.

OPERATION

- Place the massager on a flat surface, special rubber feet will prevent the device from sliding.
- Pour the water of the required temperature not exceeding the maximum level indicated on the inner side of the device. Built-in heating element is designed only to maintain water temperature, not for heating.
- Before connecting the appliance to the mains, check that the mode select switch indicates the "0" position.
- Massage roller is used for intensive massage;

DO NOT CONNECT THE DEVICE TO THE ELECTRICITY NETWORK IF YOUR FEET ARE IN THE MASSAGE BATH.

- Put your feet in the massage bath and set the desired mode using the switch:

VIBROMASSAGE (1).

- Set the mode switch to the 1 position.
- In this mode the body of the bath is vibrating allowing stimulation of the reflexogenic zones with the use of massage fingers on the bottom of the bath. This mode can be used either with water or without it. In this mode heating is not carried out.
- Vibromassage improves blood circulation and removes foot fatigue.

WARM BUBBLE BATH (3).

- Pour warm water into the massage bath and set the mode switch to the 3 position.
- In this mode compressed air passes through special holes and creates additional relaxing effect and built-in heating element maintains the water warm.
- Bubble Massage can reduce spasticity, improve the flow of oxygen into the skin of feet.

VIBROMASSAGE IN WARM WATER WITH BUBBLE BATH (2).

- For vibromassage in warm water with bubble bath set the 2 position.
- This mode is used for combined massage. Vibromassage with water heating and bubble bath improves blood circulation, relieves fatigue, stimulates nerve endings and creates relaxing effect.

DO NOT STAND ON THE MASSAGER.

CLEANING AND MAINTENANCE

- Unplug the device and let it cool completely.
- Pour water out.
- Never use abrasive detergents, organic solvents and aggressive liquids for cleaning.
- Nozzles can be washed in warm water with mild detergents.
- After massage wipe the massager with a damp cloth with mild detergent and dry it out.
- Never immerse the massager completely in water when cleaning.

DO NOT STAND ON THE UNIT.

CLEAN AND CARE

- Switch the appliance off and ensure that the massager has completely cooled.
- Do not use hard abrasives or cleaners.
- Attachments may be cleaned in soapy water
- After using, simply rinse out with water. Never immerse the entire unit into water. You can wipe off all surfaces with a cloth damped with water and a mild detergent, then rinse and dry.

STORAGE

- Be sure that the appliance is unplugged and cooled down.
- Complete instructions of section CLEAN AND CARE
- Keep the appliance in a cool, dry place.

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