SWITEL

SWISS BRAND 

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TERMOMETRO DIGITALE A CIUCCIO

BABY NIPPLE THERMOMETER

BH310

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Mode d'emploi

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Operating Instructions

D

F

I

GB

SWITEL

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Mode d'emploi 8

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Operating Instructions 18

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1 Important Safety Information

1. Carefully read the operating instructions before using the appliance.
2. Prior to taking the temperature, check that the digital thermometer, nipple and sensor are in perfect working order. Do **not** use the appliance if any damages are found.
3. The digital thermometer is solely designed for oral measuring. Temperatures to be taken under adult supervision only.
4. Do not use the digital thermometer as an ordinary nipple.
5. The digital thermometer is an electronic unit. Do not let children have it for playing.
6. Do not expose the digital thermometer to high temperatures or direct sunlight and keep it away from chemical solvents since these might get into the appliance and damage the circuitry.
7. The battery is fully sealed inside and cannot be replaced.


2 Intended Use

The digital thermometer is shaped like a nipple and is designed for orally taking the temperature of babies and infants under adult supervision. Any other use is considered unintended use. Unauthorised modification or reconstruction is not permitted. Under no circumstances open the device or complete any repair work yourself.

3 Features

- LC display for showing the temperature measured
- Memory function
- Automatic measuring and switch-off
- No mercury, no risk of poisoning by broken thermometer

4 Operating Elements

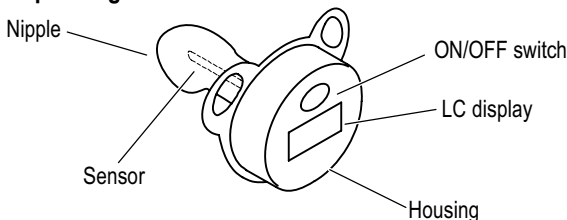5 Operation

1. To take the temperature, put the nipple into the baby's/infant's mouth.
2. Press the ON/OFF switch to turn the appliance on; the display reads 188.8 $^{\circ}\text{C}$.
The thermometer is ready for use. The temperature computed last is displayed as 37.0 $^{\circ}\text{C}$ (L0). " $^{\circ}\text{C}$ " flashes on the display.
3. If the ambient temperature is above 32 $^{\circ}\text{C}$, the temperature is shown and measuring starts.
4. If the computed temperature is above 44 $^{\circ}\text{C}$, the display reads H $^{\circ}\text{C}$.
5. If the temperature is below 32 $^{\circ}\text{C}$, the display reads L0 $^{\circ}\text{C}$.
6. When the temperature remains unchanged it is shown as $^{\circ}\text{C}$, " $^{\circ}\text{C}$ " stops flashing, and an acoustic signal is sounded ten times.
7. When measuring is complete remove the nipple from the baby's/infant's mouth.
8. Turn the appliance off to save battery power. If you do not, the appliance will turn itself off after 10 minutes. In case measuring is interrupted, start again and repeat the above steps.

6 Cleaning

1. Immerse the nipple in a 70 percent ethanol solution for 1 second.
2. Wipe the display with a cloth or use clear water.
3. The thermometer is waterproof. You can rinse it with clear water.

Attention: Do not put the thermometer into boiling water.

7 Technical Specifications

Measuring range:	32 °C - 44 °C
Measuring time:	3 minutes +/- 1 minute
Precision:	+/- 0.1 °C
Battery:	1.5 VDC
Display:	LCD (liquid crystal display)
Memory:	Displays the last temperature taken
Timer:	10x acoustic signal at stable temperature
Dimensions:	50 x 45 x 45 [mm]
Weight:	18 g with sensor sleeve
Life (battery):	approx. 100 hours

8 Disposal

In order to dispose of your device, take it to a collection point provided by your local public waste authorities (e.g. recycling centre). According to laws on the disposal of electronic and electrical devices, owners are obliged to dispose of old **electronic and electrical devices** in a separate waste container. The symbol indicates that the device must not be disposed of in normal domestic waste!

Risks to health and the environment from batteries!

Never open, damage or swallow batteries or allow them to pollute the environment. They could contain toxic and ecologically harmful heavy metals. You are legally obliged to dispose of batteries at the point-of-sale or authorised collection points where the corresponding containers are provided. Disposal is free of charge. The above symbols indicate that the batteries must not be disposed of in domestic waste and that they must be disposed of at an authorised collection point.

Packaging materials must be disposed of according to local regulations.

9 Guarantee

SWITEL equipment is produced and tested according to the latest production methods. Carefully selected materials and highly developed technology ensure perfect functioning and a long service life. The terms of guarantee do not apply to the batteries or power packs used in the products. The period of guarantee is 24 months from the date of purchase.

All deficiencies related to material or manufacturing errors within the period of guarantee will be redressed free of charge. Rights to claims under the terms of guarantee are annulled following tampering by the purchaser or third parties. Damage caused as the result of improper handling or operation, normal wear and tear, incorrect positioning or storage, improper connection or installation or Acts of God and other external influences are excluded from the terms of guarantee. In the case of complaints, we reserve the right to repair or replace defect parts or provide a replacement device. Replaced parts or devices become our property. Rights to compensation in the case of damage are excluded where there is no evidence of intent or gross negligence by the manufacturer.

If your device does show signs of a defect within the period of guarantee, please contact the sales outlet where you purchased the SWITEL device, producing the purchase receipt as evidence. All claims under the terms of guarantee in accordance with this agreement can only be asserted at the sales outlet. No claims under the terms of guarantee can be asserted after a period of two years from the date of purchase and hand-over of the product.

Declaration of Conformity

This device fulfils the requirements stipulated in the EU directive. To view the complete Declaration of Conformity, please use the download facility on our web site www.switel.com free of charge.

Hotline

In the case of technical problems, you can contact our hotline service, tel. 0900 00 1675 within Switzerland (Swisscom fees at time of going to print: CHF 2.60/min).

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This equipment complies with the European R&TTE directive.
To view the complete Declaration of Conformity, please refer to the free
download available at our web site: www.switel.com.

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In the event of technical problems, you can contact our hotline service, Tel.
0900 00 1675 within Switzerland (cost via Swisscom at time of going to print:
CHF 2.60/min).