
SWITEL

SWISS BRAND 

**DIGITALES INFRAROT-OHRTHERMOMETER
THERMOMÈTRE AURICULAIRE INFRAROUGE NUMÉRIQUE
TERMOMETRO AURICOLARE DIGITALE A INFRAROSSI
DIGITAL, INFRARED EAR THERMOMETER**

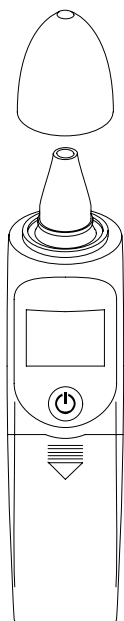BH311

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Mode d'emploi

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Operating Instructions

D

F

I

GB

SWITEL

Bedienungsanleitung	3
Mode d'emploi	9
Istruzioni per l'uso	15
Operating Instructions	21

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1 Important Safety Information

- Read the operating instructions carefully before using the device.
- Always keep the manual in the vicinity of the device and hand it over to other device users.
- Prevent exposure to smoke, dust, dirt, shocks, chemicals, moisture, heat or direct sunlight and keep it away from water, liquids and chemical solvents because these could get into the device and damage the electronics.
- Do not use the thermometer in the direct vicinity of a mobile phone. Risk of malfunction!
- Do not use the thermometer if the ear in which it is to be used has a disorder.
- Do not touch the measuring sensor with your fingers.
- The thermometer is an electronic device and must not be used as a toy by children.

2 Intended Use

The digital ear thermometer serves to measure the body temperature at the ear. Due to the short measuring time necessary and its high precision, the device is particularly suitable for measuring the body temperature of babies and small children. Any other use is considered unintended use.

Unauthorised modification or reconstruction is not permitted. Under no circumstances open the device or complete any repair work yourself.


3 Features

- Fast, precise body temperature measurement
- Last measurement is stored in the memory
- Voice output
- Automatic deactivation
- LC display
- Simple operation
- Simple cleaning
- No protective sleeve necessary for measuring sensor
- No mercury, no risk of poisoning by broken thermometer

4 Package Contents



- 1 digital, infrared ear thermometer
- 2 batteries, type AAA (LR03), 1.5 V
- 1 storage box
- 1 operating manual

5 Operating Elements

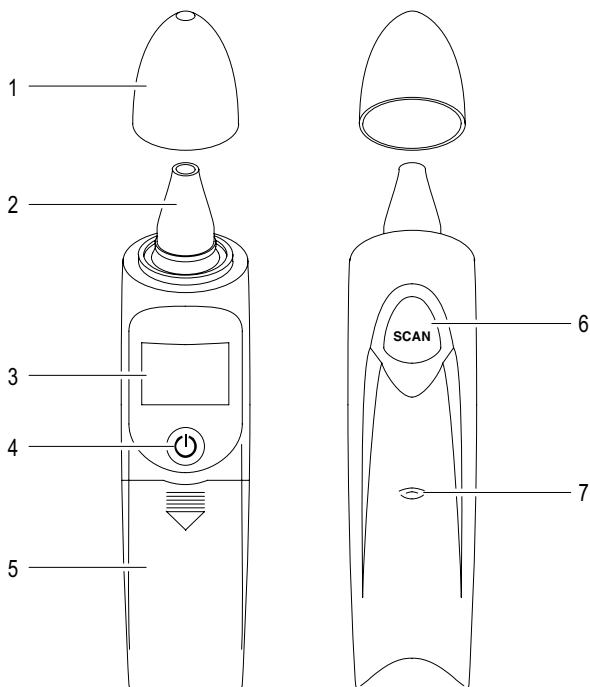- | | | | |
|---|------------------|---|---------------------------|
| 1 | Protective cap | 5 | Battery compartment cover |
| 2 | Measuring sensor | 6 | Scan button (SCAN) |
| 3 | Display | 7 | Loudspeaker |
| 4 | On/Off button | | |

Display indicators

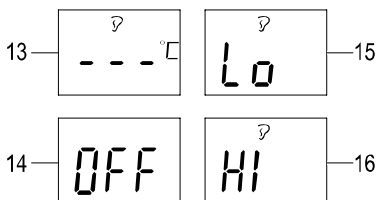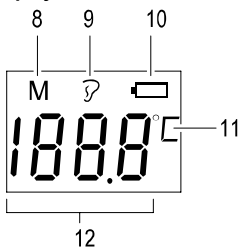
- | | | | |
|----|----------------------------------|----|---------------------------|
| 8 | Memory icon | 13 | Ready to measure |
| 9 | Ear icon | 14 | Device being switched off |
| 10 | Battery icon | 15 | Measured value < 32 °C |
| 11 | Unit of temperature measurement | 16 | Measured value > 42.9 °C |
| 12 | Display of values or information | | |

6 Inserting the Batteries

1. Slide the battery compartment cover off in the direction of the arrow.
2. Remove any batteries already inserted.
3. Insert 2 new batteries, pay attention to correct polarity, see diagram.
4. Slide the battery compartment cover back on opposite to the direction of the arrow until it audibly snaps in place.

If the top line in the display shows the battery icon and "Lo" appears, the battery charge capacity is low: change the batteries.

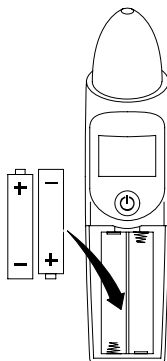You can still use the thermometer for a short time despite a low battery charge but change the batteries as soon as possible.

7 Operation

Before measuring the body temperature, ensure that

- the ear to be used has no disorders,
- the ear is free of excessive ear wax,
- the person being checked is relaxed and
- the person being checked has not just got out of a hot bath. Wait at least 30 minutes between getting out of a hot bath and measuring the body temperature.

We recommend measuring the body temperature at least 3 times and using the highest value measured as the result.

After removing the thermometer from the storage box, proceed as follows:

1. Remove the protective cap from the thermometer.
2. Press the *On/Off button* to switch the thermometer on.
 - All indicators appear briefly in the display (self-test).
 - The last temperature measured is displayed. The memory icon (M = Memory) also appears in the top line of the display.
 - The thermometer switches to Measuring mode. The ear icon in the top line of the display flashes 3 times. When the ear icon appears permanently, the thermometer is ready to use. The thermometer prompts you vocally to begin with the measurement.
3. Insert the measuring sensor into the ear:
Take hold of the ear at the top and pull the ear carefully upwards and back to stretch the naturally bent ear canal. Insert the measuring sensor carefully in the ear canal until the measuring sensor fills the ear canal.

CAUTION: Never use any force to insert the measuring sensor in the ear canal. Risk of injury!

Excessive ear wax can falsify the measured value (measured value too low). Ensure that the ear canal is free of excessive ear wax.

4. Press and hold the *SCAN button* for approx. 1 second. The thermometer measures the body temperature. Measuring is completed when two acoustic signals are sounded.
5. Remove the thermometer carefully out of the ear. Read the temperature measured in the display.
The temperature measured is also issued vocally.

If the temperature measured is under 32 °C, "Lo" appears in the display.

If the temperature measured is above 42.9 °C, "Hi" appears in the display.


If the ambient temperature is below 10 °C or above 40 °C, "Err" appears in the display.

The thermometer prompts you to complete another measurement. Repeat steps 3 to 5, even if the measurement was interrupted or is erroneous.

RECOMMENDATION: Measure the body temperature at least 3 times and use the highest value as the result.

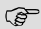6. Switch the thermometer off to save the batteries: Press the *O/Off button*.
 If the device remains unused for approx. 1 minute, it switches off automatically (automatic shutdown).
7. Clean the measuring sensor after switching off.
8. Replace the protective cap on the thermometer.

8 Cleaning

- Clean the measuring sensor with a soft cloth or cotton bud.
- Only use a soft cloth to clean the display.

9 Technical Specifications

Measuring range:	32 °C - 42.9 °C
Precision:	+/- 0.3 °C
Operating conditions:	
- Ambient temperature	Min. 10 °C, max. 40 °C
- Relative humidity	Max. 80%
Storage conditions:	
- Ambient temperature	Min. -25 °C, max. 55 °C
- Relative humidity	Max. 95%
Batteries:	2 pieces AAA 1.5 VDC (LR03)
Memory:	Display of last measured value
Dimensions:	35 x 145 x 50 [mm]
Weight:	78 g (with protective cap and batteries, without storage box)

10 Disposal

In order to dispose of your device, take it to a collection point provided by your local public waste authorities (e.g. recycling centre). According to laws on the disposal of electronic and electrical devices, owners are obliged to dispose of old **electronic and electrical devices** in a separate waste container. The symbol indicates that the device must not be disposed of in normal domestic waste!

Risks to health and the environment from batteries!

Never open, damage or swallow batteries or allow them to pollute the environment. They could contain toxic and ecologically harmful heavy metals. You are legally obliged to dispose of batteries at the point-of-sale or authorised collection points where the corresponding containers are provided. Disposal is free of charge. The above symbols indicate that the batteries must not be disposed of in domestic waste and that they must be disposed of at an authorised collection point.

Packaging materials must be disposed of according to local regulations.

11 Guarantee

SWITEL equipment is produced and tested according to the latest production methods. Carefully selected materials and highly developed technology ensure perfect functioning and a long service life. The terms of guarantee do not apply to the batteries or power packs used in the products. The period of guarantee is 24 months from the date of purchase.

All deficiencies related to material or manufacturing errors within the period of guarantee will be redressed free of charge. Rights to claims under the terms of guarantee are annulled following tampering by the purchaser or third parties. Damage caused as the result of improper handling or operation, normal wear and tear, incorrect positioning or storage, improper connection or installation or Acts of God and other external influences are excluded from the terms of guarantee. In the case of complaints, we reserve the right to repair or replace defect parts or provide a replacement device. Replaced parts or devices become our property. Rights to compensation in the case of damage are excluded where there is no evidence of intent or gross negligence by the manufacturer.

If your device does show signs of a defect within the period of guarantee, please contact the sales outlet where you purchased the SWITEL device, producing the purchase receipt as evidence. All claims under the terms of guarantee in accordance with this agreement can only be asserted at the sales outlet. No claims under the terms of guarantee can be asserted after a period of two years from the date of purchase and hand-over of the product.

Declaration of Conformity

This device fulfils the requirements stipulated in the EU directive. To view the complete Declaration of Conformity, please use the free download facility on our web site www.switel.com.

Hotline

In the case of technical problems, please contact our hotline within Switzerland under tel. 0900 00 1675 (cost with Swisscom at time of going to print: CHF 2.60/min).

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This equipment complies with the European R&TTE directive.
To view the complete Declaration of Conformity, please refer to the free download available at our web site: www.switel.com.

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In the event of technical problems, you can contact our hotline service, Tel. 0900 00 1675 within Switzerland (cost via Swisscom at time of going to print: CHF 2.60/min).